


Lab Work


CPT-based Cs microcell Atomic Clock

Marion Delehaye, Clément Lacroûte, Cyrus Rocher, Gilles Martin,
Philippe Abbé, Rodolphe Boudot



Reminder: Linear Spectroscopy (cf. Miletì, Lacroûte)

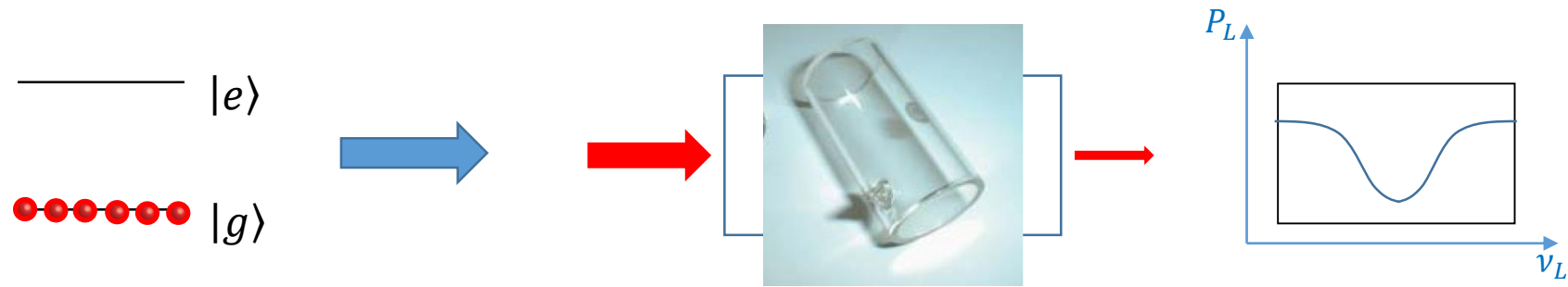

Very easy! Put a gas vapor in a glass cell!

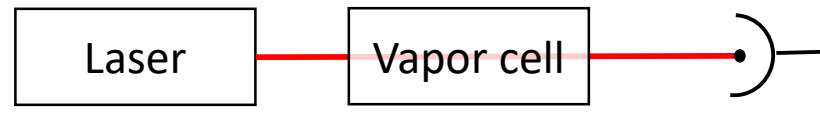




The lines are broadened by the Doppler effect

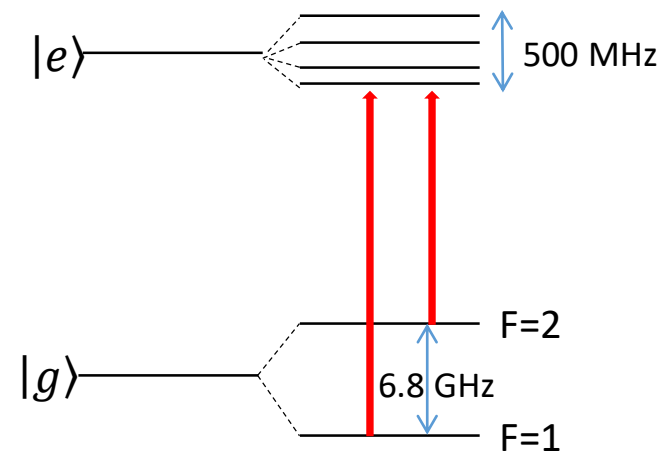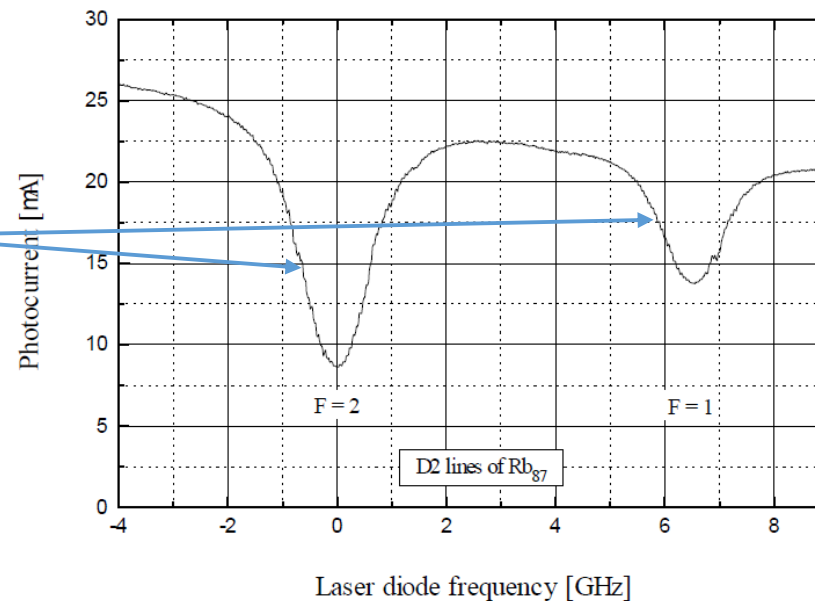

Reminder: CPT in cell clocks (cf Miletì, Affolderbach)


Coherent Population Trapping (CPT)

2 coherent light fields,
Simultaneous preparation and interrogation

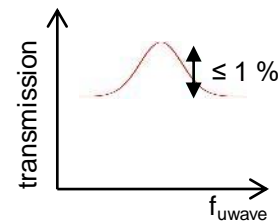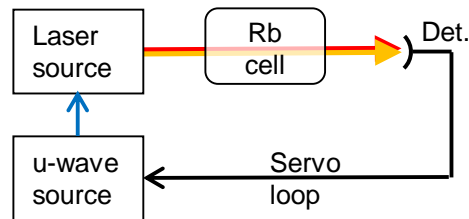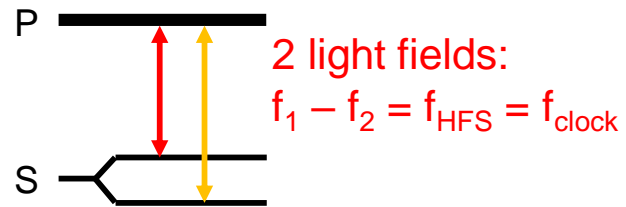

- Employs coherent laser effects
⇒ potential for new & improved clocks
(first proposal: Cyr et al., IEEE TIM, 1993)
- No microwave cavity needed
⇒ more radical clock miniaturization

Lab work CPT clock

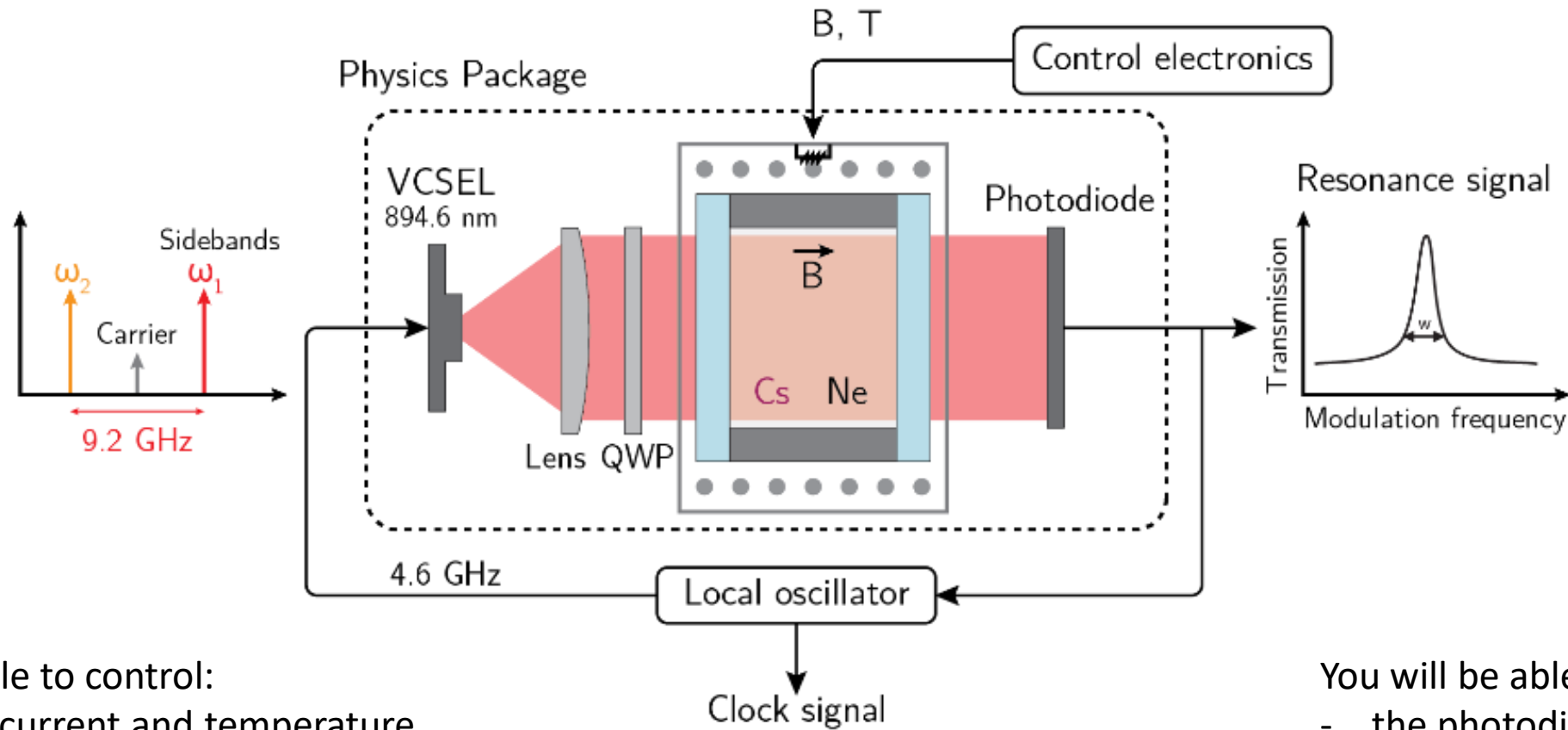

You will be able to control:


- the VCSEL current and temperature
- the magnetic field
- the local oscillator frequency
- the RF modulation power

You will be able to view:

- the photodiode signal
- laser and local oscillator error signals
- the clock signal

Lab work CPT clock – hardware


Physics package box
(VCSEL+ optics + Cs-NeMicrocell + photodiode)
Heating elements + Sensors + B-field wire+ Magnetic shielding

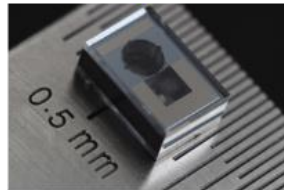

Cs-Ne microcell

Microwave coupler
(to monitor the 4.6 GHz signal)

Microwave attenuator
(manual) to tune the 4.6 GHz
power

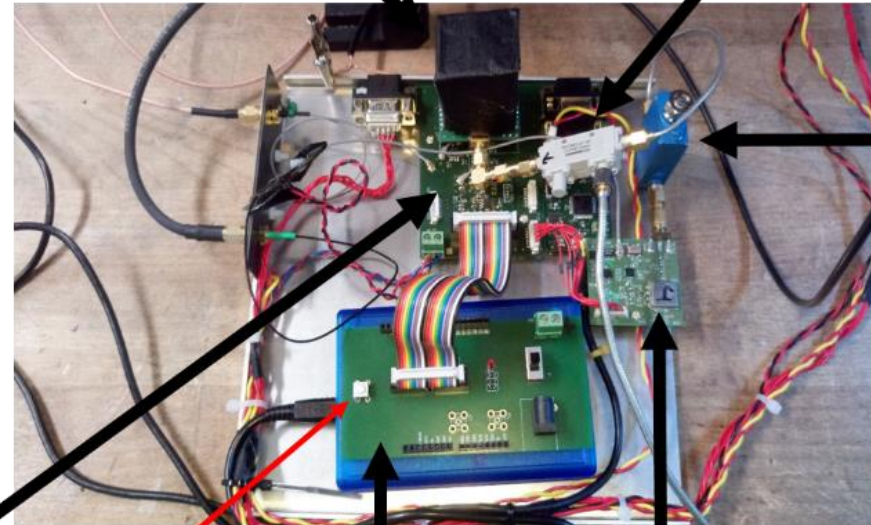


Embedded
electronics

Loops reset

Arduino

Local oscillator
(4.6 GHz VCO phase locked to
a 10 MHz quartz oscillator
using a N-fractional PLL)

Lab work CPT clock – software

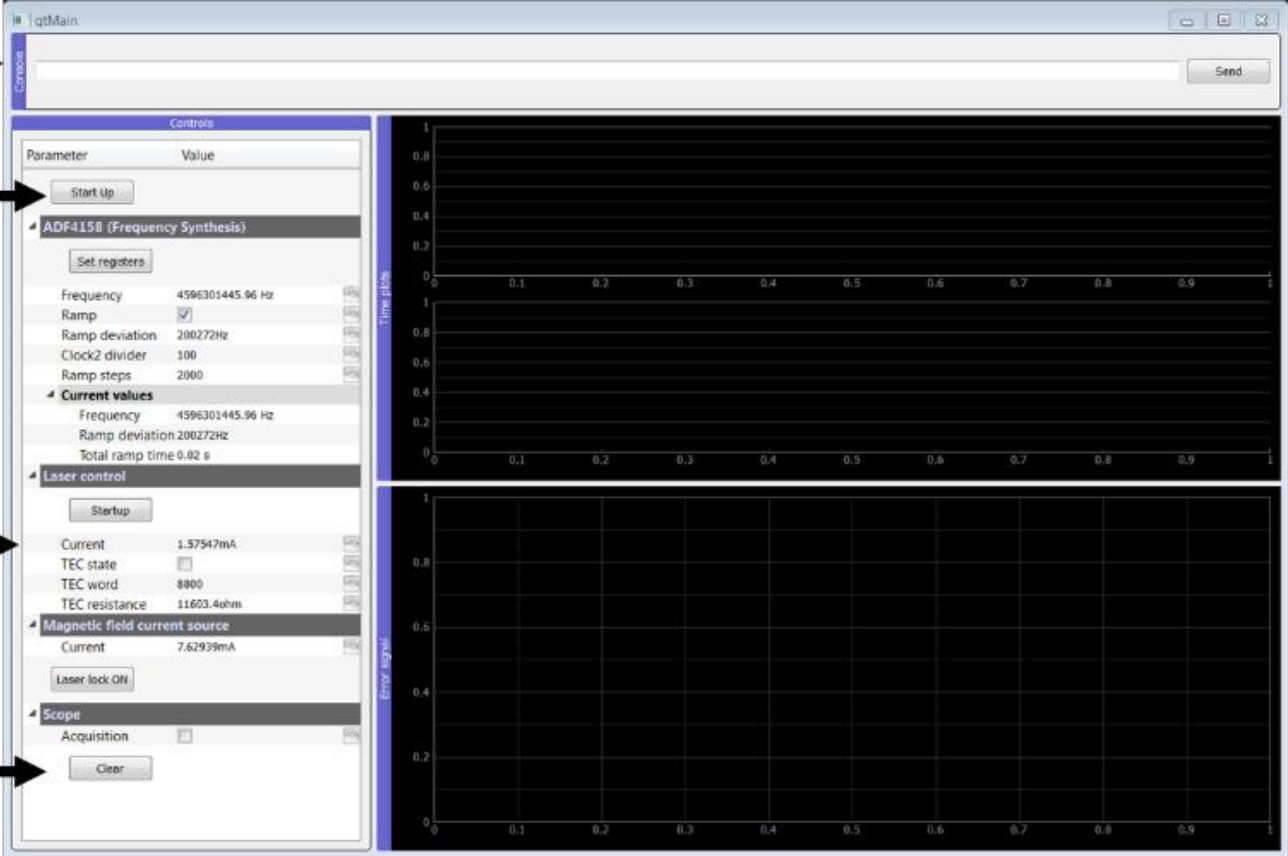

The screenshot shows the qtMain software interface. It features a top console window for serial commands, a central controls panel with various parameter settings, and three stacked plot windows on the right for real-time data visualization. Annotations with arrows point to specific UI elements: the console, the 'Start Up' button, the laser current input field, and the 'Clear' button in the scope section.

Enter serial commands here →

Startup to initialize parameters
(set the register of PLL ADFxx) →

Set the laser current here →

Clear the last scope acquisition
(use it before each new CPT scan) →

PhD output mean signal

Servo error signal

CPT error signal scan window

The whole experiment is controlled by an Arduino controller, interfaced using a C software and a Python interface. The GUI lets you input the commands as digital words within the command line. Some parameters (laser diode current, magnetic field current) can also be input in their proper unit in the dedicated boxes. You can view the commands and digital words actually sent to the micro-controller in the command window

Lab work CPT clock – what you will do

1 – Characterize the laser diode

2 – Linear spectroscopy

3 – CPT signal

4 – Lock the LO using the CPT signal

5 – Characterize the clock frequency stability

You may now start the lab work... by reading the printed labwork subject!

CPT linewidth – influence of collisions

Buffer-gas: avoids ground-state relaxation by collisions with the cell walls.

Diffusion: $\Gamma_{\text{diff}} = 2 \left(\frac{\pi^2}{L^2} + \frac{2.405^2}{R^2} \right) D_0 \frac{p_0}{p} \left(\frac{T}{T_0} \right)^{3/2} \propto \frac{1}{p}$

Collisions with the walls

Collisions: $\Gamma_{\text{coll}} = 2 n_{\text{buff}} \sigma_{\text{coll}} \bar{v}_{\text{rel}} \propto p$

Collisions with the buffer gas


Spin-exchange coll. $\Gamma_{\text{se}} = 0.6875 \bar{v}_r n_{\text{Cs}} \sigma_{\text{se}} \propto p_{\text{Cs}}$

Cs-Cs collisions

Total: $\Gamma_{12} = \Gamma_{\text{diff}} + \Gamma_{\text{coll}} + \Gamma_{\text{se}}$


\Rightarrow Compromise on $p, p_{\text{Cs}} \Rightarrow$ optimal T for a given cell (ie Cs + buffer gas mixture)

$\Rightarrow T, p, p_{\text{Cs}}$ are crucial parameters

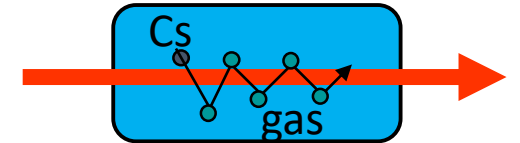

Linear absorption – linewidth


The total width of absorption lines Γ^* is determined by:

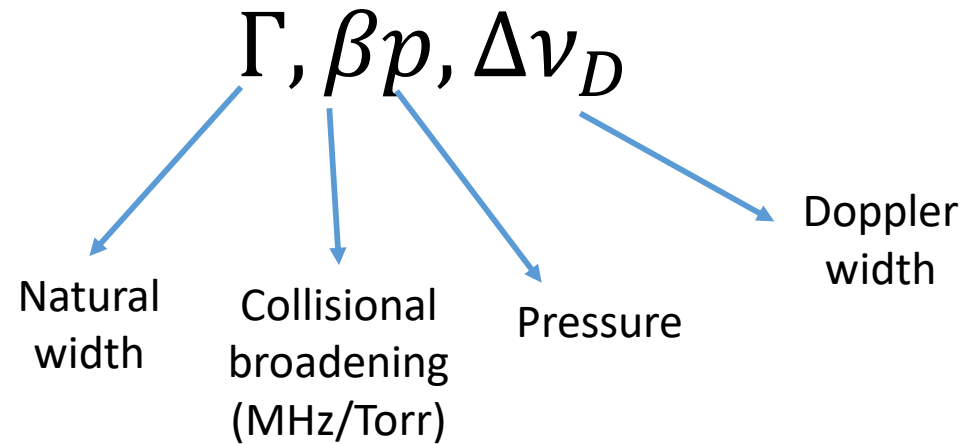

with $\Delta \nu_D = \frac{2\nu_0}{c} \sqrt{\frac{2 \ln(2) kT}{m}}$ where k is the Boltzmann constant, T is the temperature and m is the atom mass.

Linear absorption – linewidth


The total width of absorption lines Γ^* is determined by:

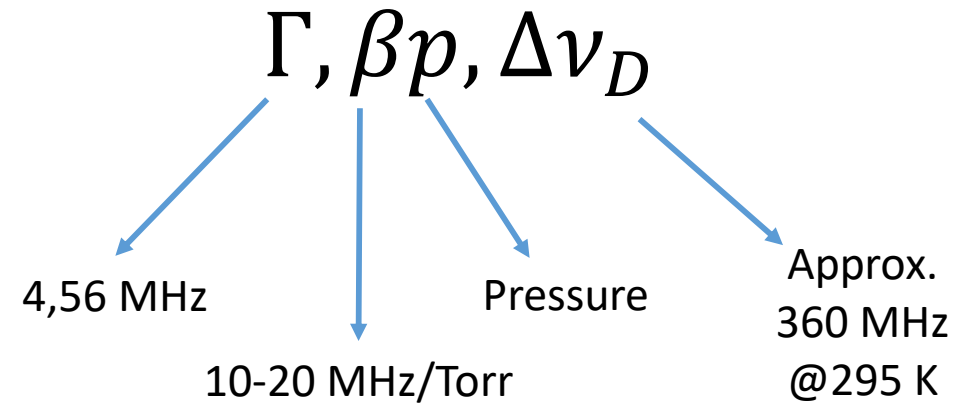

with $\Delta \nu_D = \frac{2\nu_0}{c} \sqrt{\frac{2 \ln(2) kT}{m}}$ where k is the Boltzmann constant, T is the temperature and m is the atom mass.

Linear absorption – linewidth

The total width of absorption lines Γ^* is determined by:

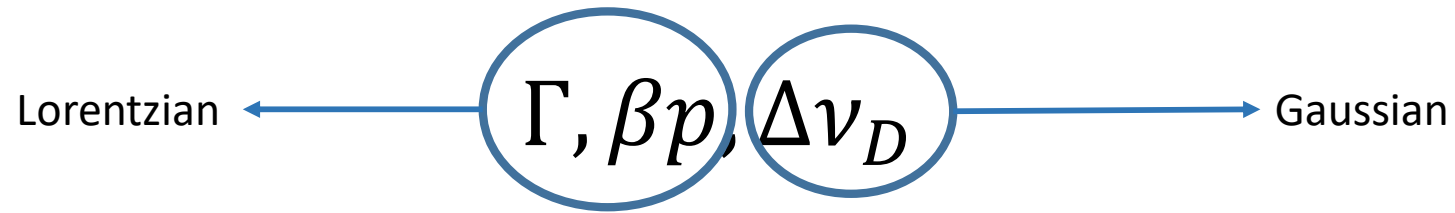

In practice, the line profile is a Voigt function (convolution of Gaussian and Lorentzian) from which one can extract both the temperature and the pressure broadening.

Linear absorption – linewidth

The total width of absorption lines Γ^* can be roughly approximated by:


$$\Gamma^* \approx \frac{\Gamma + \beta p}{2} + \sqrt{(\Gamma + \beta p)^2 / 4 + \Delta \nu_D^2}$$

From a fit, one could extract T and β or p .
Here, we know T and β so we can deduce p .

A better approximation (0.02% error – source Wikipedia):

$$f_V \approx 0.5346 f_L + \sqrt{0.2166 f_L^2 + f_G^2}.$$

Linear absorption – Beer-Lambert law

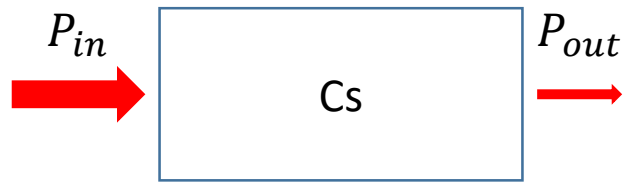

$$P_{out} = P_{in} \exp(-\alpha L)$$

where L is the cell length and α is the absorption coefficient given by

$$\alpha = \frac{\Gamma^* \lambda^3 n_{cs}}{8\pi} \sqrt{\frac{m}{2\pi kT}}$$

where n_{cs} is the Cs density inside the cell, and λ is the line wavelength.

Linear absorption – Beer-Lambert law



A diagram showing a rectangular cell labeled 'Cs'. A red arrow labeled P_{in} points into the left side of the cell, and another red arrow labeled P_{out} points out of the right side of the cell.

$$P_{out} = P_{in} \exp(-\alpha L)$$

where L is the cell length and α is the absorption coefficient given by

$$\alpha = \frac{\Gamma^* \lambda^3 n_{Cs}}{8\pi} \sqrt{\frac{m}{2\pi kT}}$$

where n_{Cs} is the Cs density inside the cell, and λ is the line wavelength.

Here we know T and Γ^* so we can deduce n_{Cs} ie p_{Cs} .

In summary:

Here we know T , L and β ;

by measuring the optical linewidth Γ^* and absorption α we can deduce the buffer-gas pressure p and the Cs pressure p_{Cs}

\Rightarrow rough estimate of expected CPT linewidth